

+ Healthcare

Sterile Services + Endoscopy

We Make Sure

Contents

Why Serchem?	3	Endoscopy	26
Where we distribute	4	Where our products fit in	28
Healthcare overview	7	Pre-treatment – bedside	30
Our specialist areas	8	Manual Sink Wash	32
Sterile Services	10	Flushing Devices	34
Where our products fit in	12	Environment Maintenance	36
Pre-treatment	14	Sink Management	37
Manual Wash	15	Accreditation	38
Ultrasonic and Automated Wash	18		
Rinse Aids	24		
Environment Maintenance	25		

+ People who rely on ultimate hygiene and decontamination can always rely on us.

Serchem is far more than just a chemical company.

Every single day we make a difference by being the most dependable, the most responsive and the most caring name in hygiene and decontamination around the world. And as a specialist manufacturer, we have complete control of everything we do, allowing us to set the ultimate standards and offer you total flexibility.

As a trusted, family business, a leading expert in Sterile Services and a market leader in Endoscopy, we have over 40 years of invaluable experience and are proud to be working with global medical distributors, private health providers as well as NHS Trusts across the UK. In that time, we've built up the most comprehensive range of more than 700 products and services to answer every challenge.

We'll draw from all of this expertise and choice and then work closely with you to understand your needs in detail and create the sharpest solutions. Whenever you need it, we can also tailor our products to meet your unique requirements with our responsive approach. Right from assessing and adapting to variations in regional water characteristics through to modifying one of our products to match your demands even more precisely.

Above all, we'll always act as a loyal partner, not simply a supplier. The unrivalled specialist who truly understands what is at stake. So, whatever we create, whatever we advise, whatever we do, **We Make Sure.**

+ Made in Britain. Distributed globally.



At Serchem we're proud of our family business values and principles which are an integral part of why our name is so trusted for delivering the very best in hygiene and decontamination. Located in the heart of England, we're equally proud to be part of the Made in Britain organisation, a growing community of businesses renowned for manufacturing to the very highest standards in Great Britain and Northern Ireland. So when you see the emblem with its forward direction arrow stamped on all of our products, it signals loud and clear the quality of everything we make and the innovative, progressive approach we take to creating the very best solutions for you.

- Chile
- China
- Germany
- Hong Kong
- India
- Ireland
- Lebanon
- Malaysia
- Malta
- Netherlands
- Oman
- Pakistan
- Philippines
- Qatar
- Saudi Arabia
- Singapore
- South Africa
- Taiwan
- Thailand
- UAE

"I have been using Serchem chemicals in various Sterile Services for more than 20 years and always found them to be value for money and very good quality."

University Hospitals of Derby and Burton NHS



+ Our service makes a world of difference.

We always stand in your shoes to understand, and solve, the challenges you face across Sterile Services and Endoscopy.

From theatres through to the sterilisation and decontamination department, Serchem is with you every step to deliver total reassurance. With premium quality, specialist solutions for any aspect of the decontamination process. With exceptional expertise. With the utmost care and attention.

And with our in-depth understanding of the complexities of endoscopes and endoscopic accessories, we've become a market leader in this highly specialised field. Which means you can always trust Serchem to clean such vitally important equipment safely using our innovative, fully compliant chemistries for the manual cleaning phase of endoscopes.

But rest assured, the highest standards don't come with a high price. Based in the UK, our significant efficiencies in production are passed on to you with very competitive costs and reliable deliveries through local distributors. And as well as working with you as a partner, we build strong, strategic alliances with respected brands to give you even greater choice and value.

So you can count on us to keep delivering. Keep developing. And keep giving you the very best service.

+ Like our chemicals, our expertise goes further.

Our experience and expertise in hygiene and decontamination for Sterile Services, Endoscopy and Dentistry allows us to offer complete solutions all across the world.

Sterile Services


Serchem is always here to take complete care of your departmental requirements. Whether you need to keep instruments moist during the transportation from theatres to the department, or you're looking for a very specific alkaline detergent safe for use on soft metals in the automated wash, you can trust us to deliver the right answers.

Endoscopy


We use our tailor-made and specifically formulated products to deal effectively with bioburden deposits found on rigid and flexible endoscopes prior to the disinfection phase. Our comprehensive range of alkaline, neutral, enzymatic and hybrid automated detergents are scientifically designed to deal with every aspect of decontamination. In particular, Serchem's enzymatic and non-enzymatic products are specifically formulated for use in Endoscopy Suites and our name is respected and trusted around the world.

Dentistry

Our extensive experience in medical hygiene and decontamination, backed by a comprehensive range of proven products, means that we can offer the same high standards of quality and service for the dental decontamination process

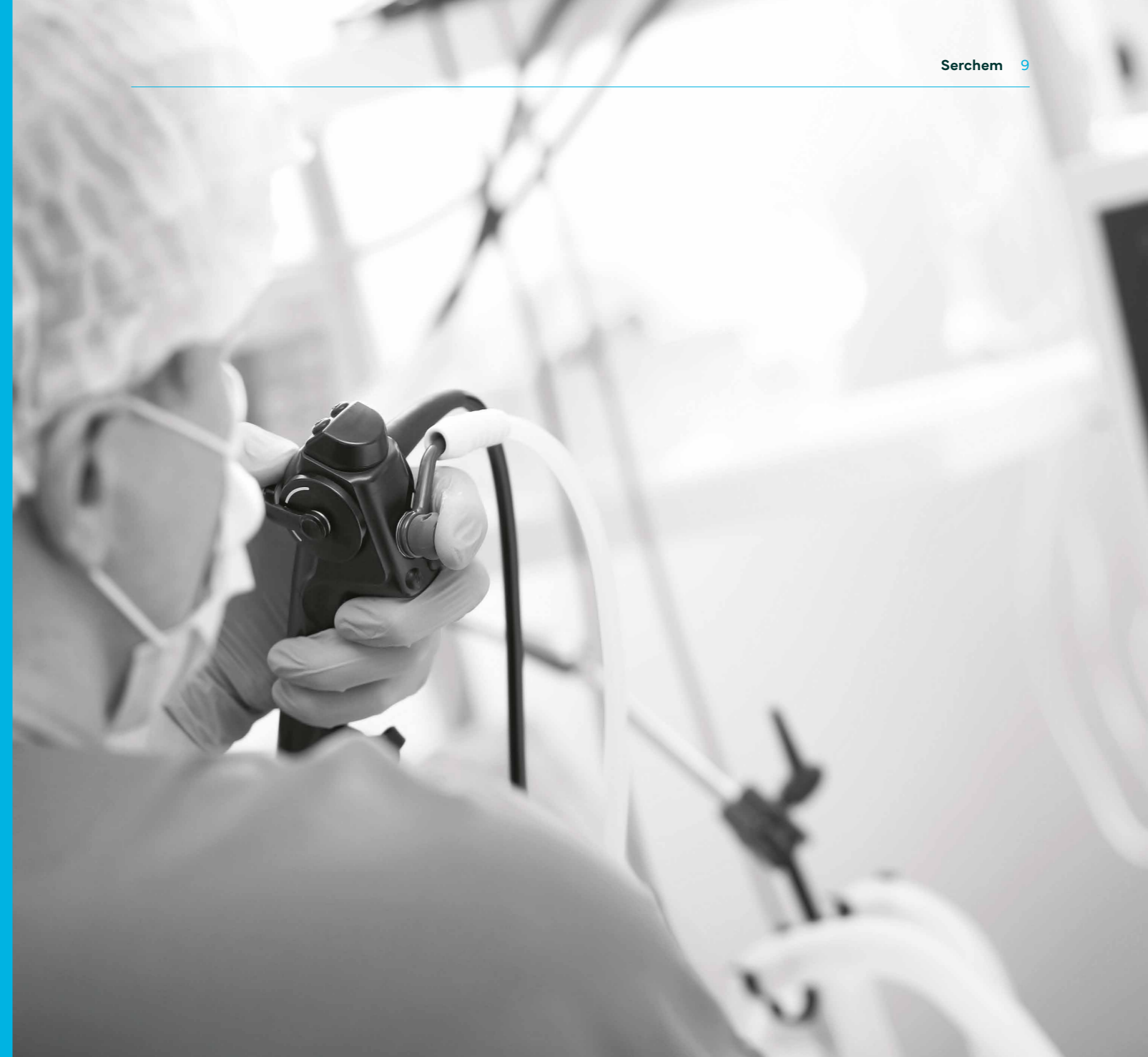

+ Sterile Services

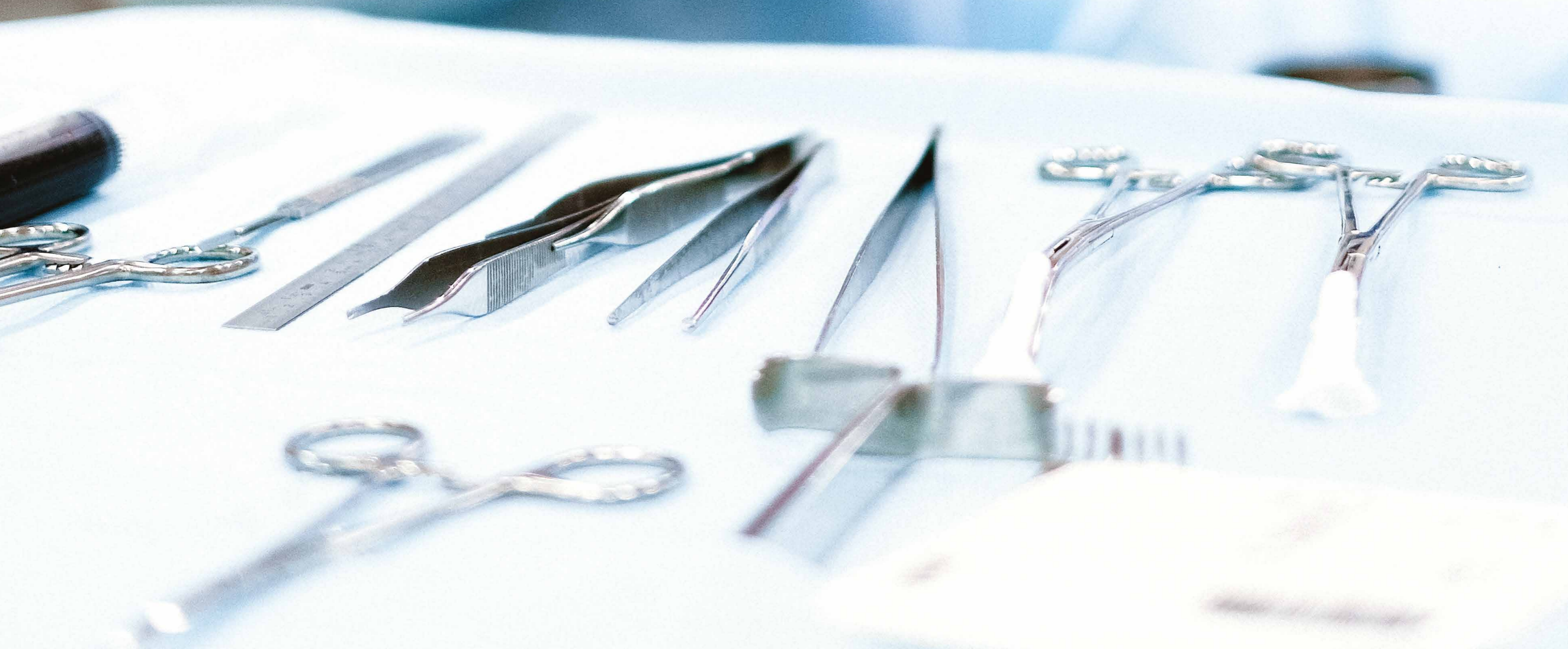

+ A perfect fit.

How Serchem fits into the decontamination sterile services process flow.



Dirty transport from theatre



Pre-treatment

Our Pre-treatment spray is designed to assist with the bioburden breakdown on surgical instruments, plastic devices and biopsy forceps.
See page 14



Manual Wash

These products are used in the sink environment to ensure that reusable medical instruments and devices such as surgical instruments are thoroughly cleaned prior to the automated process.
See pages 15–16



Ultrasonic

These products are designed for use in ultrasonic machines to assist with the removal of bioburden from surgical equipment and instrumentation.
See pages 18–22



Automatic Wash

We manufacture a range of automated detergents for use in washer disinfectors to comply with manufacturers' cleaning instructions.
See pages 18–22



Rinse Aids

We have a variety of rinse aids, both acidic and non-acidic, to support the automated process.
See page 24

Multi-Active8

 Enzymatic Formulation

Enzymatic pre-treatment foam spray for use on surgical instruments and plastic devices.


- + Multi-enzymatic pre-treatment spray.
- + Foaming nozzle on trigger spray.
- + Does not contain disinfectant.
- + Must be fully rinsed before an automated clean.



Ordering details

Pack size	6x750ml	20L
Product code	1043	104320

Triple-Zyme

 Serchem goes further
2-4ml p/ltr of wash water

 Enzymatic Formulation

Manual and machine multi-tiered enzymatic detergent for cleaning reusable medical instruments and devices.


- + pH balanced in solution.
- + Non-foaming.
- + Contains four enzymes and has three separate applications.
- + Safe for use on soft metals such as aluminium or anodised aluminium.
- + Suitable for all makes of automated washers including trolley washers.
- + Excellent cleaner in all water types.


Ordering details

Pack size	2x5L	UN2x5L	10L	20L	200L	1,000L
Product code	1088TE	1088TEUN	1088TE10	1088TE20	TE200	TE1000



PDA 51

 Serchem goes further
3-6ml p/ltr of wash water

 Alkaline Enzymatic Formulation

Alkaline hospital instrument detergent concentrate for use in washer disinfectors, ultrasonic machines, sonic irrigators and manual cleaning.

- + Alkaline.
- + Non-foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.

Ordering details

Pack size	2x5L	UN2x5L	10L	20L
Product code	PDA51/5	PDA51/5U	PDA51/10	PDA51/20



Process pH Plus

 Serchem goes further
3-6ml p/5ltr of wash water

 Alkaline Formulation

Alkaline instrument cleaner concentrate for cleaning and hand washing instruments only.

- + Alkaline sink detergent for use in the manual wash - not for endoscopes.
- + Foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Lightly fragranced.
- + Excellent cleaner in all water types.

Ordering details

Pack size	6x500ml	12x1L	2x5L	UN2x5L	10L	20L
Product code	3063TWIN	3064	3063	3063U	3068	3069



+ Always think cost per dose.



It's important to check the cost of price per dose as price per litre may not always equate to economies. Our high concentration solutions go much further and cost far less in the long run to give you real savings. This also helps you meet environmental targets with less packaging and transportation.

Serquat

General purpose cleaner concentrate for cleaning washable surfaces and hand washing instruments.

- + Neutral detergent for use in the manual wash.
- + Foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all washable surfaces.
- + Excellent cleaner in all water types.

Ordering details

Pack size	12x1L	2x5L	UN2x5L	10L
Product code	3017	3015	3015U	3016

 Serchem goes further
3ml p/ltr of wash water

 Neutral
Formulation



Manual-Active8

Super-concentrated manual enzymatic detergent for cleaning reusable medical instruments and devices.

- + pH neutral in solution, enzymatic detergent designed for use in the manual wash of surgical instruments and flexible endoscopes.
- + Non-foaming.
- + Non-fragranced.
- + Excellent cleaner in all water types.

Ordering details

Pack size	6x500ml	12x1L	2x5L	UN2x5L	10L
Product code	1039	10591	1059	1059U	105910

 Serchem goes further
4ml p/ltr of wash water

 Enzymatic
Formulation



Maxi-Active8


Super-concentrated manual enzymatic detergent for cleaning reusable medical instruments and devices.

- + pH neutral in solution, enzymatic detergent for use in the manual wash of surgical instruments, rigid endoscopes and flexible endoscopes.
- + Lightly foaming.
- + Lightly fragranced.
- + Excellent cleaner in all water types.

Ordering details

Pack size	6x500ml	2x5L	10L
Product code	1068	1052	1056


 Serchem goes further
4ml p/ltr of wash water

 Enzymatic Neutral
Formulation



"Orders are always despatched in a timely manner, expectations are realistically managed, response rate is fantastic... Serchem are an absolute pleasure to work with and we are extremely proud of the association that they have."

Distribution Partner

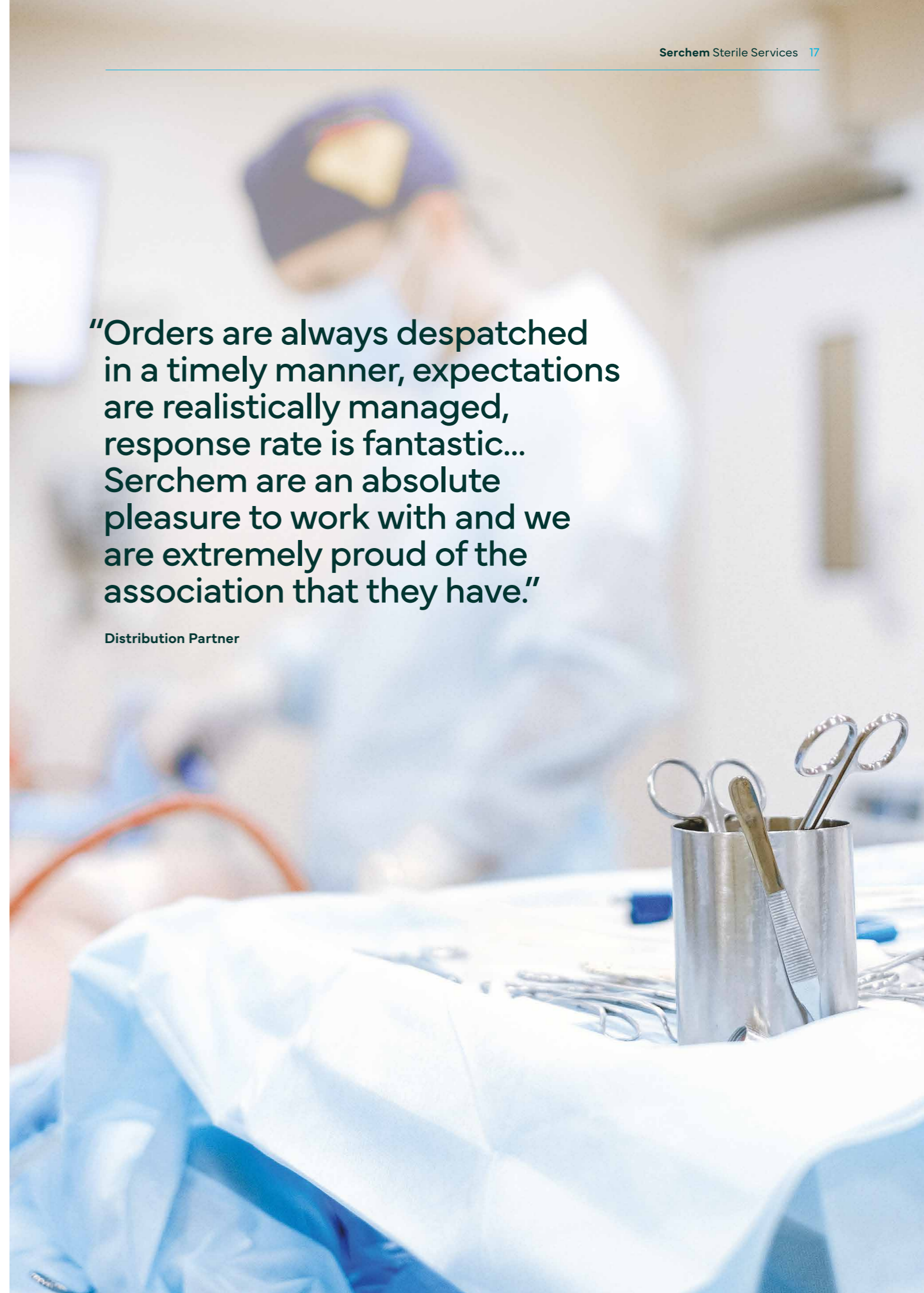

Triple-Zyme

Serchem goes further 2-4ml p/ltr of wash water Enzymatic Formulation

Manual and machine multi-tiered enzymatic detergent for cleaning reusable medical instruments and devices.

- + pH balanced in solution.
- + Non-foaming.
- + Contains four enzymes and has three separate applications.
- + Safe for use on soft metals such as aluminium or anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.



Ordering details

Pack size	2x5L	UN2x5L	10L	20L	200L	1,000L
Product code	1088TE	1088TEUN	1088TE10	1088TE20	TE200	TE1000

PDA 51

Serchem goes further 3-6ml p/ltr of wash water Alkaline Enzymatic Formulation

Alkaline hospital instrument detergent concentrate for use in washer disinfectors, ultrasonic machines, sonic irrigators and manual cleaning.

- + Alkaline.
- + Non-foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.

Ordering details

Pack size	2x5L	UN2x5L	10L	20L
Product code	PDA51/5	PDA51/5U	PDA51/10	PDA51/20



Phoenix

Serchem goes further 4ml p/ltr of wash water Enzymatic Neutral Formulation

pH balanced enzymatic detergent for use in washer disinfectors and ultrasonic machines.

- + pH balanced in solution.
- + Non-foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.

Ordering details

Pack size	12x1L	2x5L	UN2x5L	10L	20L	200L	500L	1,000L
Product code	1090	1088	1091	1092	1093	2091	51091	101091



Maximum pH Plus

Serchem goes further 2-4ml p/ltr of wash water Alkaline Formulation

Alkaline hospital instrument detergent and pre-soak concentrate for use in washer disinfectors, ultrasonic machines.

- + Alkaline.
- + Non-foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.

Ordering details

Pack size	12x1L	2x5L	UN2x5L	10L	20L	200L	500L	1,000L
Product code	1060	1078	1080	1079	1081	2078	51078	101078




+ Safe hands, from Sterile Services to Endoscopy.


As well as being a trusted name for decontamination in Sterile Services, our expertise also extends to the rigorous demands of Endoscopy cleaning. As a market leader in this highly specialised field, we use our experience, in-depth knowledge and specially formulated products to deliver complete confidence. To find out more about our Endoscopy service, please see pages 26 onwards.

Delta ES

Super-concentrated enzymatic alkaline instrument washing detergent for use in washer disinfectors.

- + Mildly alkaline.
- + Non-foaming, phosphate-free formula designed for use in the automated wash.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.

 Serchem goes further
2-4ml p/ltr of wash water

 Alkaline Enzymatic
Formulation



Ordering details


Pack size	2x5L	UN2x5L	10L	20L
Product code	107021	107021U	107020	107022




Delta

Super-concentrated alkaline instrument washing detergent for use in washer disinfectors.

- + Alkaline.
- + Non-foaming, phosphate-free formula designed for use in the automated wash.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types, especially hard water.

 Serchem goes further
2-4ml p/ltr of wash water

 Alkaline
Formulation




Ordering details


Pack size	2x5L	UN2x5L	10L	20L	200L	500L	1,000L
Product code	1070	1070U	1073	1094	2070	5070	101070

Ultra-Active8

Super concentrated machine enzymatic detergent for cleaning reusable medical instruments and devices

- + pH neutral in solution, enzymatic.
- + Non-foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of ultrasonic machines and washer disinfectors.
- + Excellent cleaner in all water types.

 Serchem goes further
4ml p/ltr of wash water

 Enzymatic
Formulation




Ordering details


Pack size	6x500ml	2x5L	UN2x5L	10L	200L
Product code	10510	1051	1051U	1053	1053200

Auto-Active8

Super-concentrated machine enzymatic detergent for cleaning reusable medical instruments and devices.

- + pH neutral in solution, enzymatic.
- + Non-foaming.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of ultrasonic machines and washer disinfectors.
- + Excellent cleaner in all water types.

 Serchem goes further
4ml p/ltr of wash water



 Enzymatic
Formulation



Ordering details

Pack size	6x500ml	2x5L	UN2x5L	10L	20L	200L
Product code	1069	1054	1054U	1055	1058	1054200

Progenica

 Serchem goes further
2-4ml p/ltr of wash water  Alkaline
Formulation

Ultra-concentrated alkaline washing detergent for use in washer disinfectors and ultrasonic machines.



- + Alkaline.
- + Non-foaming, phosphate-free formula designed for use in the automated wash.
- + Safe for use on soft metals such as aluminium and anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.



Ordering details

Pack size	12x1L	2x5L	UN2x5L	10L	20L	200L	500L	1,000L
Product code	1098	1095	1095U	1096	1097	2095	5095	101095

Bravo TW2

 Serchem goes further
2-4ml p/ltr of wash water  Acidic
Formulation

Acidic trolley washing detergent for use in the automated trolley wash.

- + Specifically designed to deal with the issue of soft anodised transport trolleys.
- + Non-foaming.
- + Non-fragranced.
- + Balances the wash water in decontamination process to reduce pH.
- + Aqueous solution dispensed in usual way.



Ordering details



Pack size	10L	20L
Product code	108910	1089

“Serchem provides us with brilliant and responsive service at all times. Their products are high quality with reliable and repeatable results. They have a solution (literally) for all cleaning requirements or processes and accommodating water hardness/chemistry.”

Nuffield Health HSSU



Alpha CA

 Serchem goes further 1-2ml p/5ltr of rinse water  Acidic Formulation

Citric acid-based rinse agent and neutraliser for laboratory ware and instrument washer disinfectors.



- + Acidic rinse aid designed for use in all types of automated cleaning of hospital instruments.
- + Neutralises any alkaline residues.
- + Suitable for stainless steel instrumentation.
- + Does not cause corrosion or oxidation.
- + Suitable for all water conditions.



Ordering details

Pack size	2x5L	UN2x5L	10L
Product code	1071CA	1071CAU	1077CA

Sequel

 Serchem goes further 1-2ml p/5ltr of rinse water  Neutral Formulation

Acid-free rinse agent for instrument washer disinfectors and trolley washers.

- + Low foaming.
- + For use on all surgical instruments and can be used on trolley washers.
- + Non-fragranced.
- + Used in the final rinse stage of the automated decontamination process.
- + Aids drying process of reusable surgical instrumentation.



Ordering details

Pack size	12x1L	2x5L	UN2x5L	10L	20L	200L	500L
Product code	10741	1074	1075	1076	1085	2074	2075

Clean Zone

Biocidal sanitiser spray for use on hard surfaces.



- + Lightly foaming.
- + For use on tables, trolleys and walls and all surfaces that come into contact with instrumentation.
- + Very effective hard surface cleaning.
- + Particularly effective against MRSA.
- + Economical.



Ordering details

Pack size	6x750ml
Product code	1064

Renovate Plus

 Serchem goes further 5-10ml p/ltr of wash water  Acidic Formulation

Acidic-based instrument enhancing detergent for removing rust and corrosion from stainless steel surgical instruments.



- + Lightly foaming.
- + Designed for stainless steel surgical instruments. Not for use on anodised/aluminium instruments.
- + Can be used in the cleaning/decontamination cycle to clean washer disinfectant and internal tubes to prevent residual and water scales from sticking to the chamber walls.



Ordering details

Pack size	12x1L	2x5L	UN2x5L
Product code	1066	1065	1065U

Scale Away

 Serchem goes further 100ml p/5ltr of wash water  Acidic Formulation

Descaling compound for the removal of limescale.

- + Powerful removal of limescale
- + Safe to use on stainless steel equipment.
- + Economical.



Ordering details

Pack size	12x1L	2x5L
Product code	1061	4003

+ Together we always deliver.

Reliability is in our nature. We have to be completely dependable as a team and across absolutely everything we do. It's the core of our business. From our products to our people to the service we give you, we deliver consistency and quality, day in, day out. And you can always rest assured that when you need an even more responsive answer, it's always done with the utmost care, attention and rigour.



+ Endoscopy

+ A trusted part of the process.

How Serchem fits into the decontamination endoscopy process flow.



Endoscopy cleaning



Pre-treatment – bedside

We have a range of convenient non-enzymatic bedside packs and sponges specifically designed for the cleaning of flexible and rigid endoscopes. These can be supplied in several different formats, depending on requirements. **See pages 30–31**



Manual Sink Wash

These products are used in the sink environment to ensure that flexible and rigid endoscopes are thoroughly cleaned prior to the automated process. **See pages 32–33**



Flushing Devices

These products simplify pre-cleaning, improve efficiency and promote compliance to regulatory standards as well as help with operator fatigue. The dosing equipment also reduces human error and the risk of overusing chemicals. **See pages 34–35**

Flexi-Clenz Waterless Kit Eco



Non-enzymatic bedside pack in compostable pots for manual cleaning of flexible and rigid endoscopes.

- + Easy tear sachets ensure safer cleaning process for patient and end user.
- + Lightly foaming.
- + Suitable for all makes of flexible and rigid endoscopes.
- + Lightly fragranced.
- + Suitable for all water conditions.
- + Eco friendly compostable pot.

Ordering details

Pack size	1x250
Product code	1030.3



Flexi-Clenz Waterless Kit Plastic



Non-enzymatic bedside pack in plastic pots for manual cleaning of flexible and rigid endoscopes.

- + Easy tear sachets ensure safer cleaning process for patient and end user.
- + Lightly foaming.
- + Suitable for all makes of flexible and rigid endoscopes.
- + Lightly fragranced.
- + Suitable for all water conditions.

Ordering details

Pack size	1x250
Product code	1030.2



Flexi-Clenz Sachets



Non-enzymatic detergent in easy tear sachet for manual cleaning of flexible and rigid endoscopes.

- + Easy tear sachets ensure safer cleaning process for patient and end user.
- + Lightly foaming.
- + Suitable for all makes of flexible and rigid endoscopes.
- + Lightly fragranced.
- + Suitable for all water conditions.

Ordering details

Pack size	1x250
Product code	1036CS.2



+ Always think cost per dose.



It's important to check the cost of price per dose as price per litre may not always equate to economies. Our high concentration solutions go much further and cost far less in the long run to give you real savings. This also helps you meet environmental targets with less packaging and transportation.

Cylindrical Sponges

Sponges with a cylinder-shaped design for easy grip.

- + Durable, high quality medical sponge.
- + Available to order as a separate item.
- + Can be used to protect the endoscope distal end during transit and in the pre-cleaning phase.
- + Dimensions – OD 50mm x ID 10mm x L 80mm

Ordering details

Pack size	1x480
Product code	1029



Neutra-Scope

 Serchem goes further
4ml p/ltr of wash water

 Neutral
Formulation

**Neutral manual and machine detergent
for cleaning reusable medical instruments and devices.**


- + Non-foaming.
- + Great cleaning power.
- + Suitable for all makes of flexible and rigid endoscopes.
- + Odourless.
- + Suitable in all water conditions.

Ordering details

Pack size	2x5L	UN2x5L	10L	20L
Product code	1034	1034U	1035	103420



Flexi-Clenz

 Serchem goes further
2-4ml p/ltr of wash water

 Neutral
Formulation

**Neutral manual and machine detergent
for cleaning reusable medical instruments and devices.**


- + Non-enzymatic neutral detergent for manual cleaning of flexible and rigid endoscopes.
- + Lightly foaming.
- + Suitable for all makes of flexible and rigid endoscopes.
- + Lightly fragranced.
- + Suitable for all water conditions.

Ordering details

Pack size	2x5L
Product code	1036C



Neutra-Scope AF

 Serchem goes further
4ml p/ltr of wash water

 Neutral
Formulation

**Neutral manual and machine detergent
for cleaning reusable medical instruments and devices.**

- + Alcohol free.
- + Non-foaming.
- + Suitable for all makes of flexible and rigid endoscopes.
- + Odourless.
- + Suitable in all water conditions.

Ordering details

Pack size	2x5L	UN2x5L
Product code	1034AF	1034AFU



Manual-Active8

 Serchem goes further
4ml p/ltr of wash water

 Enzymatic
Formulation

**Super-concentrated manual enzymatic detergent
for cleaning reusable medical instruments and devices.**


- + pH neutral in solution, enzymatic detergent designed for use in the manual wash of surgical instruments and flexible endoscopes.
- + Non-foaming.
- + Non-fragranced.
- + Excellent cleaner in all water types.


Ordering details

Pack size	6x500ml	12x1L	2x5L	UN2x5L	10L
Product code	1039	10591	1059	1059U	105910



Triple-Zyme

 Serchem goes further
2-4ml p/ltr of wash water

 Enzymatic
Formulation

**Manual and machine multi-tiered enzymatic detergent
for cleaning reusable medical instruments and devices.**

- + pH balanced in solution.
- + Non-foaming.
- + Contains four enzymes and has three separate applications.
- + Safe for use on soft metals such as aluminium or anodised aluminium.
- + Suitable for all makes of automated washers.
- + Excellent cleaner in all water types.


Ordering details

Pack size	2x5L	UN2x5L	10L	20L	200L	500L
Product code	1088TE	1088TEUN	1088TE10	1088TE20	TE200	TE10001



Maxi-Active8

 Serchem goes further
4ml p/ltr of wash water

 Enzymatic Neutral
Formulation

**Super-concentrated manual enzymatic detergent
for cleaning reusable medical instruments and devices.**

- + pH neutral in solution, enzymatic detergent for use in the manual wash of surgical instruments, rigid scope and flexible endoscopes.
- + Lightly foaming.
- + Lightly fragranced.
- + Excellent cleaner in all water types.

Ordering details

Pack size	6x500ml	2x5L	UN2x5L
Product code	1068	1052	1056



ScopeTech™ DF

The ScopeTech™ DF 4-in-1 flexible endoscope cleaning system combines a host of benefits. This highly versatile unit delivers everything from scope inflation for accurate leak testing to precise detergent delivery, two flushing/ rinsing points as well as scope inflation and cleaning process reports. This step-by-step cleaning process tracks every scope that's cleaned. The system also promotes compliance with reprocessing requirements while reducing the time it takes to reprocess a scope.



Ordering details

Product code SCOPEDF

FlushTech™ GI

The FlushTech™ GI eliminates the need for manual syringing. This dual port, high capacity flushing system promotes compliance with endoscope flushing requirements while reducing the time it takes to reprocess a scope. This is an easy-to-operate, robust system for flushing detergent, rinse water and air through all internal channels and any elevator wire channels of a flexible endoscope. The system also features reporting capabilities.



Ordering details

Product code 7260030

EndoFlate™ GI

EndoFlate™ GI is an automated inflation system used to pressurise the flexible endoscopes for leak testing. It provides an easy-to-read digital screen with a numerical pressure value, timer, pressure bar and date/time. A red screen indicator notifies the technician if the liner fails to meet desired inflation pressure for an immersion leak test. In the event of a minor leak, EndoFlate™ GI can keep the scope inflated so the technician can clean the scope without risk of fluid ingress.



Ordering details

Product code 7260040

Scope-Clenz

Flushing solution specifically formulated for use in the FlushTech™ GI/ScopeTech™ DF flushing devices only.

- + Low hazard sanitising solution.
- + Clear in solution.
- + Odourless.
- + Safe in application for all members of staff.
- + Supplied in the ready to use 300ml or 5L container.



Ordering details

Pack size	2x5L	30x300ml Pouch	30x300ml Bottle
Product code	1001	1002P	1002

Clean Zone

Biocidal sanitiser spray for use on hard surfaces.

- + Lightly foaming.
- + For use on tables, trolleys and walls and all surfaces that come into contact with instrumentation.
- + Very effective hard surface cleaning.
- + Particularly effective against MRSA.
- + Economical.




Ordering details


Pack size	6x750ml
Product code	1064

Scale Away

Descaling compound for the removal of limescale.

- + Powerful removal of limescale
- + Safe to use on stainless steel equipment.
- + Economical.

 Serchem goes further
100ml p/5ltr of wash water

 Acidic
Formulation



Ordering details

Pack size	12x1L	2x5L
Product code	1061	4003

“I cannot put into words the support that the Serchem team have provided. I have dealt with Serchem for about 10 years now and can honestly say from service to advice they are outstanding.”

Distribution Partner

+ Sink Management.

These Sink Management solutions for Sterile Services and Endoscopy have been developed to deliver excellent economical usage with their ease of application.

Working together, the carefully thought through elements offer a highly efficient and common-sense approach to assist technicians actively in the manual cleaning phase of the decontamination process and ensure a safe working environment. In addition, all of the following equipment can be both installed and maintained by Serchem.

DoseTech™ DMS

The DoseTech™ DMS eliminates all the variability of bottle mounted pumps by simply delivering precise amounts of detergent into the sink. It includes built-in programmable settings, allowing users to set dosing parameters that will meet specific conditions for removing soil at various sink levels for achieving optimum cleaning results.



Product code	7165068 – DoseTech DMS with no USB or flow meter	Product code	7165069 – DoseTech DMS with flow meter
Product code	7165070 – DoseTech DMS with flow meter and USB		

Q10+ Pump System

With its intuitive and easy-to-use interface, the Q10+ single peristaltic pump system delivers a consistent measure of detergent or rinse aid for every cleaning process. In addition, configuration and set up are both quick and simple.

Product code	Q10
--------------	-----



Temperature Sensor

With its high visibility digital display, this easy to use and essential piece of equipment makes it completely simple for technicians to achieve and maintain the correct water temperature during the manual cleaning phase of the decontamination process.

Product code	1109
--------------	------



Sink Level Indicator

We'll work closely with your department to ensure the required water level is precisely achieved. We'll then install the Sink Level Indicator to assist technicians in working to the correct detergent/water dilutions to achieve optimum manual cleaning results. We can also calibrate the DoseTech DMS and the Q10+ units to the department's specified water capacity to ensure economical 'in-use' costings.

Product code	1108
--------------	------



+ A dependable name backed by respected accreditation.

When it comes to quality, expertise and true value for money, you can always be sure of Serchem. And this pledge to you is backed by accreditation from a wide range of highly respected organisations.



+ We Make Sure.

To discover the Serchem difference and how much further we go to give you the ultimate in hygiene and decontamination, make sure you get in touch today. We would welcome the opportunity to talk about your requirements and how we could help.

Serchem

Heatherset, Arleston, Telford, Shropshire, TF1 2LY, UK

T +44 (0)1952 223130

E support@serchem.co.uk

serchem.com

